

Independent suture strands distribute tension on the repair to create a larger footprint

The high-strength solution for shoulder repairs...

- **46% stronger** than CorkScrew*, which pulls out at 59 lbs. (vs. 86 lbs. for FASTIN RC)¹
- **NOW pre-loaded with ORTHOCORD, ETHIBOND® Excel or PANACRYL® suture**
- **Excellent soft bone purchase** achieved with smaller inner core and wide cancellous threads
- **Reduce suture-binding** with independent suture eyelets that allow for a modified Mason-Allen technique
- **Easy single-step insertion** with sharp tip in either single or double-row repairs

Now pre-loaded with
ORTHOCORD Suture
with needles

FASTIN[®] RC

Dual-Channeled Anchor

FASTIN[®] RC

Dual-Channeled Anchor



Laser lines and letters for easy identification of suture color and eyelet orientation during insertion



Drill point requires less force to insert than the tapered tip design


Anchor's wide threads and small core optimize bone purchase

Unique T-Hex design provides independent suture channels to help prevent suture binding

Two sutures dispense the load over more of the tendon

Average Tensile Strength¹

FASTIN RC Anchor vs. Corkscrew* Anchor

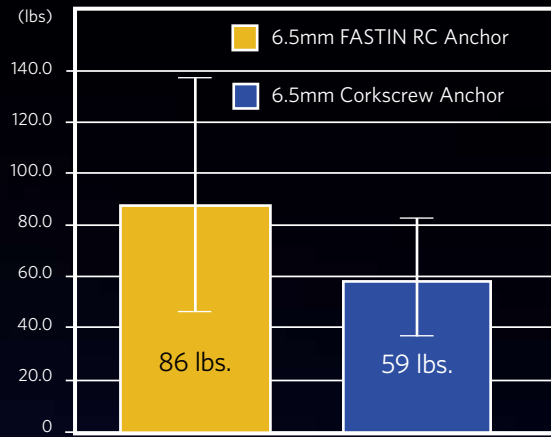

Error bars indicate +/- one standard deviation from the mean

ORDERING INFORMATION

NUMBER	DESCRIPTION	QUANTITY	ORDER QTY
222993	FASTIN RC 5mm ORTHOCORD™ w/o Needles		
222723	FASTIN RC 5mm ETHIBOND® w/o Needles		
222743	FASTIN RC 5mm PANACRYL® w/o Needles		
222994	FASTIN RC 6.5mm ORTHOCORD w/o Needles		
222764	FASTIN RC 6.5mm ETHIBOND w/o Needles		
222765	FASTIN RC 6.5mm PANACRYL w/o Needles		
222991	FASTIN RC 5mm ORTHOCORD w/ Needles		
222722	FASTIN RC 5mm ETHIBOND w/ Needles		
222742	FASTIN RC 5mm PANACRYL w/ Needles		
222992	FASTIN RC 6.5mm ORTHOCORD w/ Needles		
222762	FASTIN RC 6.5mm ETHIBOND w/ Needles		
222763	FASTIN RC 6.5mm PANACRYL w/ Needles		

Please consult package insert for indications, contraindications, precautions, warnings and adverse effects.



Dyneema® Purity is a registered trademark of Royal DSM N.V.

For more information, call your DePuy Mitek representative at 1-800-382-4682 or visit us at www.depuymitek.com. DePuy Mitek, Inc., 325 Paramount Drive, Raynham, MA 02767

© DePuy Mitek, Inc., 2007. All rights reserved. Printed in the USA. P/N 900951 6/07

* Corkscrew is a trademark of Arthrex Inc. ¹ Testing done in cadaveric specimens. Data on file.

